# HDC-300C Titan HD Camera Instructions for Use Version 1.0

Sunoptic Technologies®  
6018 Bowdendale Avenue  
Jacksonville, FL 32216 USA

Customer Service: 904 737 7611  
Toll Free 877 677 2832

**EC|REP** AJW Technology Consulting GmbH  
Königsallee 106  
40215 Düsseldorf, Germany  
[+49 \(0\) 211 3013 2232](tel:+4921130132232)

## TABLE OF CONTENTS

1. INTRODUCTION
  - 1.1 INDICATIONS FOR USE
  - 1.2 FUNCTION OF DESIGN
  - 1.3 SPECIFICATIONS
  - 1.4 CLASSIFICATIONS
2. GENERAL WARNINGS
3. PRODUCT OVERVIEW
  - 3.1 THE CAMERA CONSOLE
  - 3.2 THE CAMERA HEAD
  - 3.3 COMPATIBLE ACCESSORIES
4. SETUP AND INTERCONNECTION
  - 4.1 SETTING UP THE CONSOLE
  - 4.2 CONNECTING THE CAMERA
  - 4.3 CONNECTING MONITORS
  - 4.4 CONNECTING AC POWER
  - 4.5 POWERING ON THE UNIT
  - 4.6 TURNING OFF THE UNIT
5. OPERATION
  - 5.1 USING THE CONSOLE TOUCH SCREEN INTERFACE
  - 5.2 ENTER PATIENT DATA
  - 5.3 REVIEW/PLAY CASES
  - 5.4 RECORD
  - 5.5 SETTINGS
  - 5.6 MAIN MENU
  - 5.7 KEYBOARD CONTROLS
  - 5.8 FOOTSWITCH
  - 5.9 FAULT DISPLAYS
6. CLEANING AND DISINFECTION
7. MAINTENANCE, SERVICING & REPAIR, WARRANTY
  - 7.1 USER SERVICE AND MAINTENANCE
  - 7.2 LIMITED WARRANTY
  - 7.3 REPAIR
8. END OF PRODUCT LIFE
9. SYMBOLOGY

## 1. INTRODUCTION

Congratulations on the purchase of your new HD Camera System!

The Sunoptic Technologies High Definition Camera System is a head-mounted surgical camera used to capture still and video images of surgical procedures.

You have chosen the best and we would like to make sure you receive optimal results with your new Camera System by using it correctly.

This Operator Manual will help you to install, operate and clean the Camera System. It will give you maintenance and service guidelines as well as recommendations for best performance results.

Federal law (United States of America) restricts this device to use by, or on the order of, a physician.

### 1.1. INDICATIONS FOR USE

The intended use of this device is to display, record and play back surgical procedures. As such, this device is not intended to diagnose, prevent, monitor, treat or alleviate disease, injuries or handicaps, nor for the investigation, replacement or modification of the anatomy or of a physiological process, and is not to be used with or as an endoscope or related system. As such, this device has no essential performance in accordance with IEC 60601 and users are advised that should this system fail to operate, suitable backup systems be available for any procedure in which this device may be used.

There are no known contraindications.

The HD Camera is intended to be used in a controlled operating room environment with compatible devices by qualified medical personnel. The camera system and head are provided non-sterile.

### 1.2. FUNCTION OF DESIGN

The HD Camera System incorporates a remote camera head which captures an image as presented with a coupled light source and displays this image onto a viewing monitor. Displayed images and videos may be captured and stored internally or transferred or transmitted via a variety of means, controlled through the device's integrated touch panel.

### 1.3. SPECIFICATIONS

PARAMETER	VALUE	
Power Requirements	Voltage: 90-240VAC Frequency: 47-31Hz Power: 100VAC/2.9A 240VAC/1.27A	
Video Outputs	2x HDMI 2x DVI	
Connectivity	2 x USB 2.0 ports (The USB port on the rear outer edge is used only for software updates) 2 x Ethernet ports	2 x Camera Inputs (Rear Camera Input Optional) 1 x Footswitch Input
Internal Storage (Typical)	12Mbps = 90 hours 18Mbps = 60 hours 24Mbps = 45 hours	
CCU Dimensions	13" W x 15.25" L x 5.25" T	
CCU Weight	16.0 lbs	
Camera Weight	0.85lbs	
Transport & Storage Conditions	-40 to 122 <sup>o</sup> F (-40 to 50 <sup>o</sup> C) 10 to 90%, non-condensing	
Operating Conditions	50 to 86 <sup>o</sup> F (10 to 30 <sup>o</sup> C) 30 o 75% RH, non-condensing	

### 1.4 CLASSIFICATIONS

PARAMETER	VALUE
System Classification	FDA Class I EU Class I Health Canada Class I
Design Criteria	IEC 60601-1 Ed. 3.1: 2012
EMC Certifications	CISPR 11 Class B IEC 60601-1-2 Ed. 4
CE Marking	CE Marking for MDD 93/42/EEC
According to protection against electric shock.	Class I [grounded]
According to Degree of protection against harmful ingress of water	Camera Control Units are Ordinary [IPX-0] no protection.
According to the degree of safety in the presence of Flammable Anesthetics	Equipment is NOT suitable for use in the presence of flammable anesthetics.
According to the mode of operation	Continuous

## 2. GENERAL WARNINGS

Use of this equipment may present hazards to the user and/or end patient. Before operating this device, please read this operating manual thoroughly and follow all warnings, cautions, and instructions for use. The words warning, caution, and note carry special meaning and should be carefully reviewed:

**Warning:** Indicates risks to the safety of the patient or user. Failure to follow warnings may result in injury to the patient or user.

**Caution:** Indicates risks of improper use and/or damage to the equipment. Failure to follow cautions may result in loss of function or product damage.

**Note:** Indicates special information to clarify instructions or present additional useful information.

 An exclamation mark within a triangle is intended to alert the user to the presence of important operating and maintenance instructions in the manual.

### **Warnings**

 To avoid potential serious injury to the user, patient, and/or damage to this device, please note the following warnings:

1. Failure to follow instructions in this manual may lead to serious injury or damage to the equipment. Read this operating manual thoroughly, especially the warnings, and be familiar with its contents before using this equipment.
2. This equipment is designed to be used by a qualified physician, having complete knowledge of the use of this equipment and the procedure to be performed.
3. This equipment should be installed and tested prior to use. A pre-operative check should be performed prior to administration of patient anesthesia to ensure all desired functions are operational, that a viable surgical image is being displayed on the monitor, and there are no signs of damage to the equipment.
4. It is recommended to have a back-up system ready for use in the event of primary equipment failure.
5. This equipment may present a risk of electric shock. To reduce this risk, this equipment must only be connected to a supply mains with protective earth.
6. In the case of equipment failure, this equipment may draw excessive power from the supply circuit and interrupt service to other equipment powered by the same circuit. To reduce this risk, this equipment should not share the same electrical outlet or grounding with life support or life sustaining equipment.
7. An uninterruptable power supply is recommended.
8. Always set up the console in a location that allows adequate ventilation (airflow) to the console. Insufficient ventilation may cause the console to overheat and shut down and create risk of fire. Always set up the console to position the appliance inlet (where the power cord attaches to the equipment) so that it is readily accessible.
9. Use only the provided medical grade power cord and replacement fuses specified by the manufacturer to reduce risk of harm resulting from fire or interruption of mains circuit. Follow the instructions provided herein and always remove power from the device prior to inspection or replacement of fuses.

10. Use only manufacturer's specified compatible accessories and peripheral equipment. Use of non-approved equipment may cause loss of function.
11. Before each use, check the outer surface of this equipment to be used to ensure that there are no rough surfaces, sharp edges, or protrusions that can cause injury.
12. Avoid dropping the camera system or rough handling. The camera system contains sensitive components prone to damage by mechanical shock.
13. This device complies with the IEC 60601-1 safety standard. When peripherals are connected to this device, a Medical Electrical (ME) System is formed and the system should be evaluated for conformance to IEC 60601-1 safety standards. When used with other equipment, the leakage currents may be additive. The person creating the ME system is responsible to comply with the applicable safety regulations and safety standards for their location. Connect only IEC 60601-1 approved equipment to this device. Never touch peripheral equipment connections on this device and the patient at the same time. This can create a risk of shock to the patient.
14. Do not use the equipment in the presence of flammable liquids, gases, or other materials susceptible to ignition due to electrical sparking.
15. This equipment generates and radiates RF energy which may affect the normal function of nearby installed equipment. Equipment that generates and radiates RF energy may affect the normal function of the HD Camera System. When choosing a location for the HD Camera System, consult the "Electromagnetic Compatibility" section of this manual to ensure proper function with other installed equipment.

## **Cautions**

To prevent improper use and/or damage to this device, please note the following cautions:

1. Carefully unpack this unit and check if any damage occurred during shipment. If damage is detected, refer to the Warranty and Return Policy section of this manual.
2. This equipment generates heat and uses cooling fans during normal use. To reduce the risk of overheating and possible loss of function and/or damage to the equipment, install the equipment in a well-ventilated space within the operating conditions listed herein.
3. There is a risk of damage to the equipment if the equipment is subjected to cleaning or sterilization methods not approved by the manufacturer. To reduce risk of loss of function and/or damage to the equipment, only use the approved cleaning and sterilization methods described herein.
4. High frequency surgical devices may cause interference with the video image displayed on the surgical monitor. To reduce this risk, high frequency devices and imaging equipment should be located on separate supply circuits.
5. There are no user serviceable parts inside the console or camera heads. There are hazardous voltages present inside the console. Do not remove the cover. Return the device to the manufacturer for service.

**NOTE:** The warranty is void if any of these warnings or cautions are disregarded

### 3. PRODUCT OVERVIEW

#### 3.1 THE CAMERA CONSOLE

Figure 1 - Camera Console Front Panel

No.	Name	Function
1	Touch Screen	Controls all aspects of Camera Console
2	ON/OFF Button	Turns the unit on or off
3	USB	Removable storage or keyboard connectivity
4	Camera Port	Accepts the camera cable

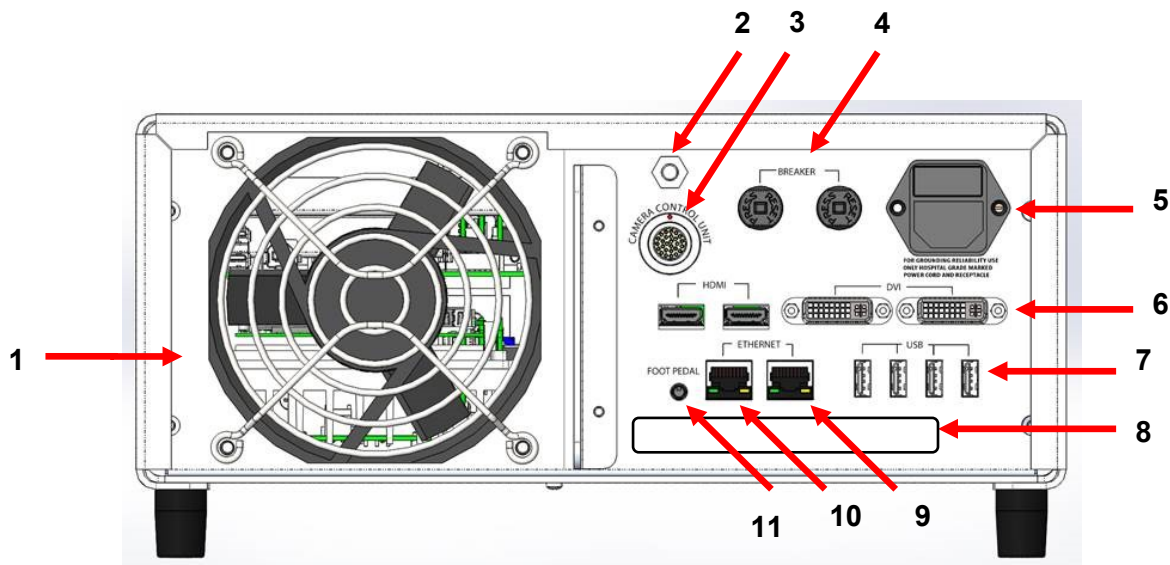
Figure 2 - Camera Console Rear Panel

No.	Name	Function
1	Fan	Provides air flow for cooling of unit
2	Grounding Connector	For electric potential equalization
3	Camera Port	Accepts the optional second camera cable
4	Circuit Breakers	Over current protection
5	AC Power Entry and Switch	Accepts AC Power cord and switches on AC power
6	Video Output Ports	Two HDMI provide stabilized video and two DVI ports provide unprocessed video
7	USB Input Ports	Port functions from left to right 1: Keyboard or removable storage connectivity 2: No Function 3 and 4: Used for firmware updates, not to be used for keyboard or file transfer.
8	Label	Regulatory Product Label
9	Recorder Ethernet Port	Allows user to connect to recorder via Ethernet
10	Stabilizer Ethernet Port	Allows user to connect to stabilizer via Ethernet
11	Foot Pedal Port	Accepts the 2.5mm plug from foot pedal

### 3.2 THE CAMERA HEAD

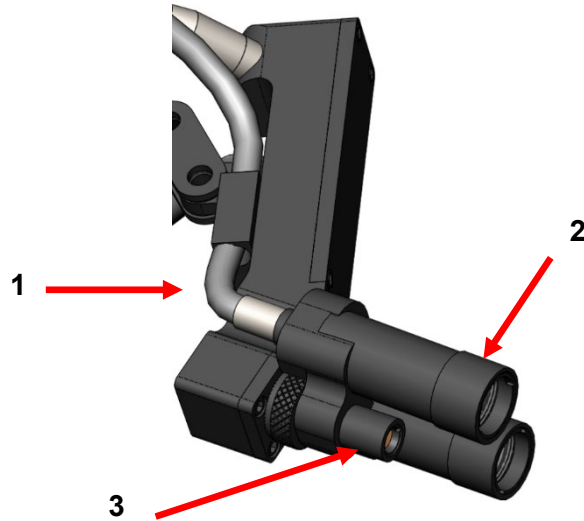
Figure 3 - Camera Head

No.	Name	Function
1	Fiber Optic Cable	Transmit light from external light source
2	Headlight Module	Projects light onto surgical site
3	Camera Lens	Lens through which video is recorded

### 3.3 COMPATIBLE ACCESSORIES

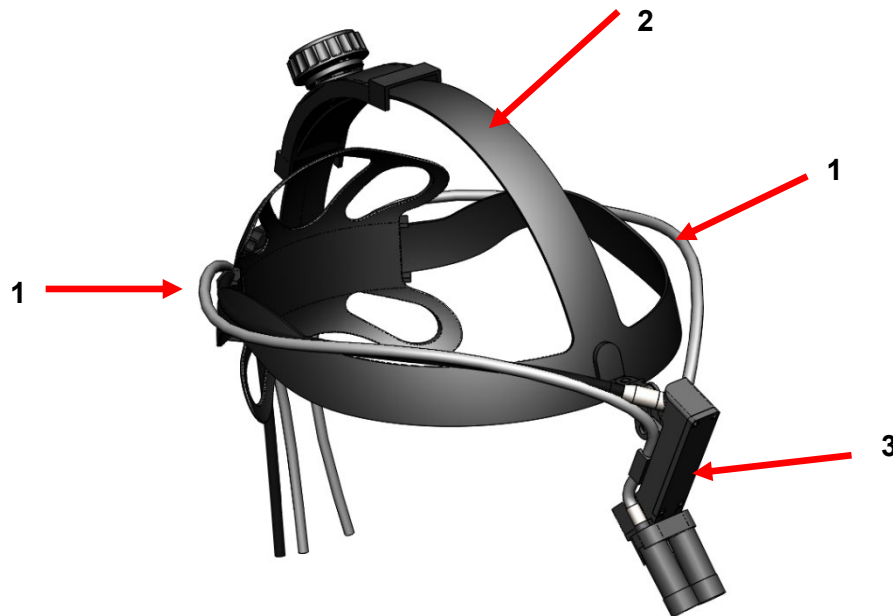
Figure 4 - Compatible Accessories for the Camera System

No.	Name	Function
1	Fiber Optic Cables	Transmit light from external light source
2	Headband	Adjustable headband system
3	Camera Head	Record video

## 4. SETUP AND INTERCONNECTION

### 4.1. SETTING UP THE CONSOLE

Place the Camera Console on a stable surface (cart, counter, stand etc.).

NOTE: Avoid places where the console may be splashed with liquid. Absolutely DO NOT use in any environment with explosive or flammable gases. This camera console pulls cooling air into its cabinet through a grille in the bottom panel underneath the front portion and exhausts warm air from the fan in the rear.

DO NOT block the space in front or below the front of the camera console or impede exhaust air flow behind or above the cabinet.

A Clear space minimum of 5"/12.7cm behind and above the cabinet is required. The console should not be placed where its exhaust will influence other devices, nor where exhaust from other devices will influence the unit.

### 4.2. CONNECTING THE CAMERA

**WARNING:** THE CAMERA CORD HAS A LOCKING CONNECTOR. TO REMOVE THE CABLE FROM THE UNIT, GRAB THE CONNECTOR AND PULL THE CONNECTOR OUT. THE CONNECTOR CANNOT BE REMOVED BY PULLING THE CABLE AND DOING SO MAY RESULT IN A DAMAGED CABLE (THIS TYPE OF DAMAGE IS NOT COVERED BY THE PRODUCT WARRANTY).

Ensure that the camera cable is properly inserted into the Camera Port. To do this, rotate the right-angle connector before inserting until the red dot on the connector aligns with the red dot on the receptacle. Push in firmly.

### 4.3. CONNECTING MONITORS

The Sunoptic Technologies High Definition Camera System comes equipped with four (4) video outputs (two (2) HDMI and two (2) DVI). All four (4) outputs can be used at once and are configured as shown below.

Output	Configuration
HDMI 1	Stabilized Video
HDMI 2	Stabilized Video
DVI 1	Unstabilized Video
DVI 2	Unstabilized Video

Follow the instructions for setting up the monitor(s) and connect to the Sunoptic HDC-300 Camera Console using either HDMI or DVI cables.

#### 4.4. CONNECTING AC POWER

Ensure the rear power switch is in OFF position.

Connect the AC power cord to the power inlet located on the rear panel of the Camera Console. Ensure the power cord is completely seated into the power inlet and engage the cord retainer if so equipped.

Use only power cords provided with the HD Camera System.

Plug the AC power cord into a wall outlet using the three (3) prong plug supplied with the unit.

To prevent electric shock, connect power cords of peripheral equipment through medical isolation transformers.

**Note:** When using medical isolation transformers, be sure to check the transformer power ratings. Ensure that the power cord is connected to the main power with three prong plugs (USA use UL2601-1 rated isolation transformers and/or power strips only).

**CAUTION:** Nearby radio or TV units may subject this equipment to radio interference. To avoid adverse electromagnetic effects, DO NOT operate this equipment near any RF energy equipment.

#### 4.5. POWERING ON THE UNIT

After ensuring that all external devices are properly connected, use the switch on the power entry module on the back of the unit to turn the power ON and the switch will light up GREEN. Once AC power is applied to the unit, push the On Button on the front panel to turn the unit ON. The Sunoptic splash screen will appear on the touch screen after turning the unit on.

Please allow approximately 30 seconds for the system to boot.

#### 4.6. TURNING OFF THE UNIT

After pressing the button on the front of the unit to turn the Camera Console off, allow up to a full minute for the system to shut down before switching the unit off at the power entry switch or trying to unplug the unit.


Please allow approximately 30 seconds for the system to shut down completely.

## 5. OPERATION


### 5.1 USING THE CONSOLE TOUCH SCREEN INTERFACE

Once the unit has been powered on and finishes its boot sequence the **Home Screen** will appear. The following options are available.

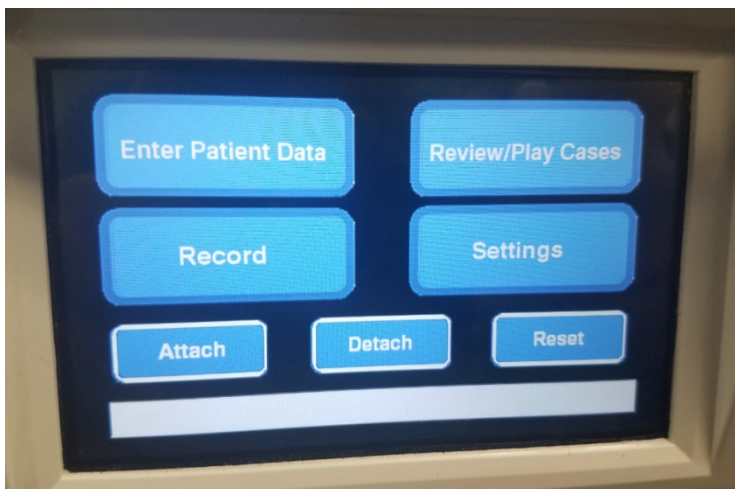
### 5.2 ENTER PATIENT DATA

This menu option will detail instructions on how to access the Patient Data entry form within the camera system. Add information using the external keyboard. Press Enter to return to begin recording video or Esc to stop data entry.

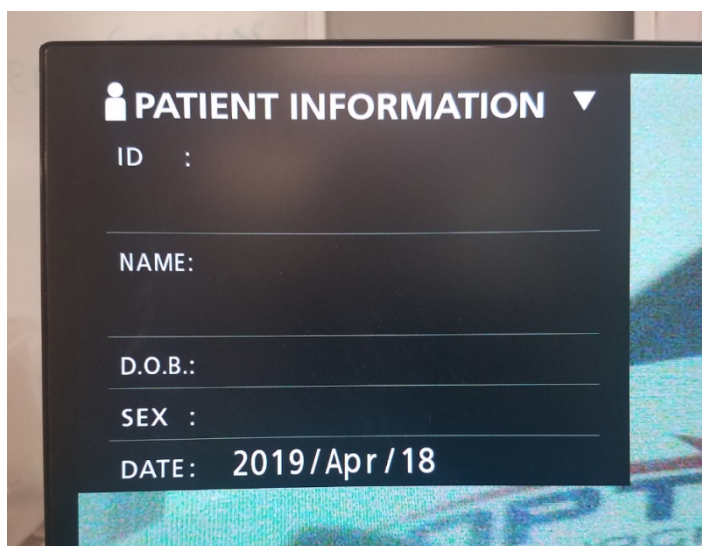
### 5.3 REVIEW/PLAY CASES

Use this menu to watch previously recorded cases, with included play, pause, fast forward and stop options while in video playback. "Go to Prev File" and "Go to Next File" takes the user to the previous or the next video or picture within the selected case. Press "Close file and exit" to stop video playback and return to the main screen.
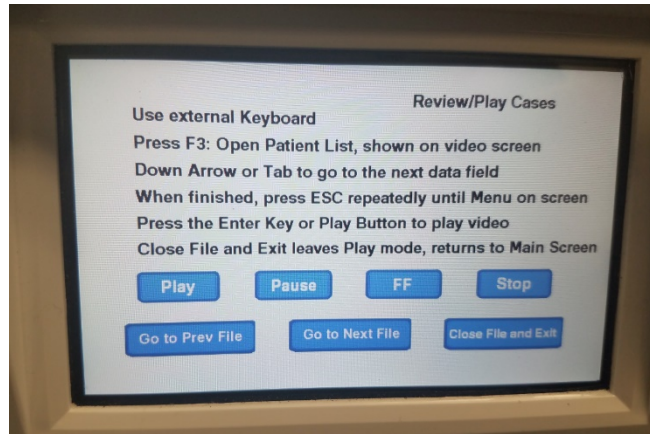
### 5.4 RECORD

The Record Button gives the operator access to recording features of the camera system including record, pause, and transferring files to a USB drive.

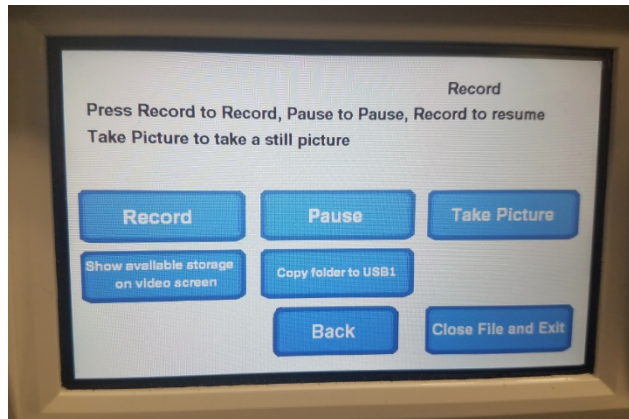Record	Starts video recording (if patient data is entered, Case will be stored with name, otherwise Case will be automatically named with the date and time video starts).
Pause	Pauses video recording (pressing this button a second time will not start recording again), will start at this point when commenced by pressing record.
Take Picture	Takes a static image of whatever is on the screen at the time.
Show available storage	Toggles the option to show (in hours and minutes) the remaining internal memory and the free space of the USB drives inserted into the front ports.
Copy Folder to USB1	Copies the selected file to the USB drive inserted into USB Port 1.
Back	Exits Record menu to main menu, video recording will continue.
Close file and Exit	Exits Recording file to Home Screen, video recording will stop.

## 5.5 SETTINGS

Access this menu to see the Product number or current Firmware version in the unit. Contact Sunoptic Technologies to use the Update button in this menu.

## 5.6 MAIN MENU

Press F4 on the keyboard from the Home Screen to access the Main Menu for the Camera System.

Menu Option	Choices	Descriptions
Settings	Current Settings	Shows all current settings for the camera system.
	Playback From	Selects source for video playback (choose from iHDD or USB1).
	Date and Time	Change the date and time manually by using the arrows on the keyboard.
	Language	Change the language on the recording device.
Photo	Quality	Select desired quality for pictures.
	Format	Select desired format for pictures.
Video	Quality	Select desired quality for videos. - EQ (Economy Quality) = lowest quality picture, uses least memory - SQ (Standard Quality) = average quality picture - HQ (High Quality) = highest quality, uses most memory
	Auto Index	Select the length of time that indices are set at for playback within the fast forward function.
Advanced	Patient Info. Entry	Switch between entering patient information with a keyboard or by using a barcode scanner. <b>DO NOT SELECT 'OFF' OR ALL FUNCTIONS OF THE KEYBOARD WILL BE LOST AND THE UNIT WILL HAVE TO BE REPROGRAMMED AT THE FACTORY.</b>

	Auto File Division	Set the size for automatic file division at either 4 GB or 1 GB. 1 GB is standard to make file transfer easier.
	Format Recording Device	Format the recording device (FORMATTING A DEVICE WILL ERASE ALL RECORDED DATA, BE AWARE THAT THIS CANNOT BE UNDONE).
	OSD	Set how information is shown on an external monitor. ON = all information is shown ICON = only status icon is shown ICON + TIMER = only timer and status icon is shown OFF = information is not shown
	Copy to (from Internal HDD)	Copy recorded videos/photos from the internal Hard Drive to a different location. Choose between starting from the newest case folder or the oldest. USB1: Copies information to USB 1 Drive Network (FTP): Contact Sunoptic Technologies for additional information. Network (CIFS): Contact Sunoptic Technologies for additional information.
	Erase iHDD (Keep)	Choose what data is erased from the iHDD, data will be erased when the unit is started. Folders will be erased starting with the oldest data to keep either 200GB, 100GB, 50GB, the last 24 hours of data or the last 7 days of data. Select the OFF option to avoid automatically deleting any data, this will not allow recording once the iHDD is full.
	Storage Device Info.	This shows various information about the selected storage device. (S.M.A.R.T. options for each device will provide a memory diagnostic for the device)
	Beep	This feature will beep whenever the internal recorder loses power and hasn't been told to turn off, leave this feature on at all times.
	Printer Options	This function is not available currently.
	Firmware Update	Sunoptic Technologies will issue a notice and a file on a flash drive with instructions for any future firmware updates.
Network	Mode Select	See "Copy to (from Internal HDD)". Contact Sunoptic Technologies for additional information.
	IP Address	Choose between automatic and manual IP Address.
	Current IP Address	Shows all current IP Address information.
Exit	-	Exits Main Menu.

## 5.7 KEYBOARD CONTROLS

The keyboard included with the camera system can be connected to the unit via the front panel USB port or the inner most USB port located on the rear panel of the console.

Key	Function
F3	<p>Brings up case list, scroll with arrow keys.</p> <p>Press F4 to bring up secondary menu within case list, giving options of Copy, Edit, Delete and Cancel. Press enter on the option required. Copy will open the copy menu, choose the desired destination for the file and press enter (file transfer will take some time to complete.) Edit will allow the user to edit the patient information for the case that is selected. Fill in the new information and press enter to make the changes to the case name. Select delete to remove the file from the internal memory and Cancel to exit this menu. Pressing F4 again will exit this menu. Once out of the secondary menu, press ESC or F3 again to exit the original case menu.</p> <p>Press enter on a case name to access all files associated with it. For a file with multiple indices OR videos and pictures, use the up and down arrows to access the additional files.</p> <p>Press F12 to access Record &gt; Close File and Exit to exit playback.</p>
F4	<p>Brings up Main Menu for the camera system (see section 5.6). Press ESC or F4 to exit this menu.</p>
F5	<p>Press F5 to immediately begin recording a session. Press F5 again to pause. Press F6 to stop video, pressing F5 again will start recording in the same case under a new index. Press F12 or access Record &gt; Close File and Exit to stop recording.</p>
F6	<p>Stops the index that is being recorded during an active video.</p>
F7	<p>Takes a still picture.</p>
F8	<p>Enters Playback system. Use the up and down arrows to access the additional files within a selected case. Press F12 or access Record &gt; Close File and Exit to stop playback and return to home screen.</p>
F10	<p>Enters Patient Data menu to edit the patient information for a new case. Fill in the new information and press enter to start recording the first index. Press F12 or access Record &gt; Close File and Exit to stop recording.</p>
F11	<p>Same functionality as F3.</p>
F12	<p>Closes file and exits recording for any video playing or being recorded.</p>

## 5.8 FOOTSWITCH (OPTIONAL EQUIPMENT)

Use the footswitch in conjunction with the unit to control video recording and playback and to take pictures without accessing the keyboard or the touch screen. Plug the footswitch cable into the 2.5mm jack on the back of the unit.

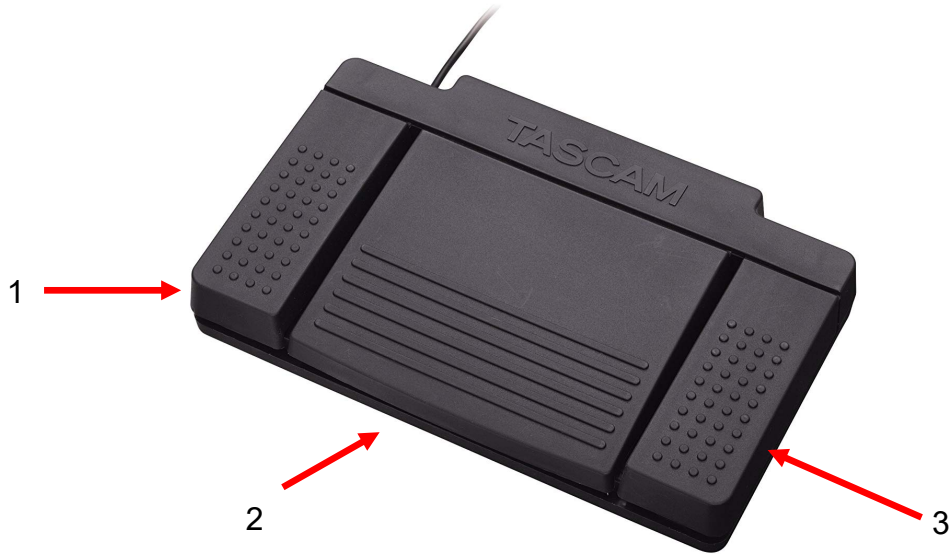
Figure 5 - TASCAM Footswitch

No.	Name	Function
1	Stop	Stops video recording or playback.
2	Video/Pause	Starts new video or new file, pauses video that is playing or recording.
3	Photo	Take photo of current screen.


## 5.9 FAULT DISPLAYS

There are two fault conditions detectable by sensors that will render the Camera System inoperable until they are corrected. A fault condition notice will appear on the screen identifying the fault and how to correct it.

### FAN FAULT


In the event of a fan problem the unit will display this screen.

Turn off and restart the Camera System. If the fault persists the unit needs servicing.
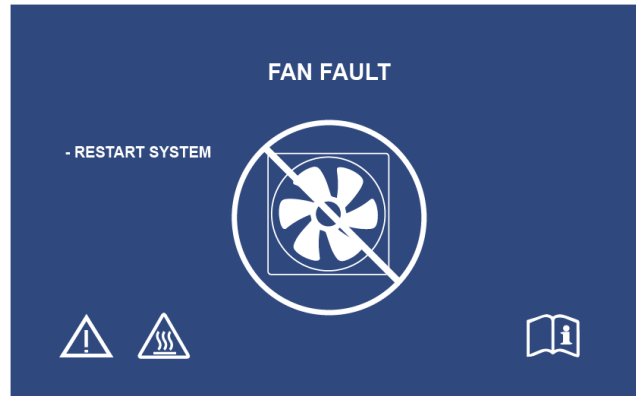
### SYSTEM IS COOLING

If the internal air temperature exceeds preset limits this screen will display. The fan will run to reduce the temperature. If acceptable reduction is achieved within ten seconds, the system will turn back on to resume normal operation.
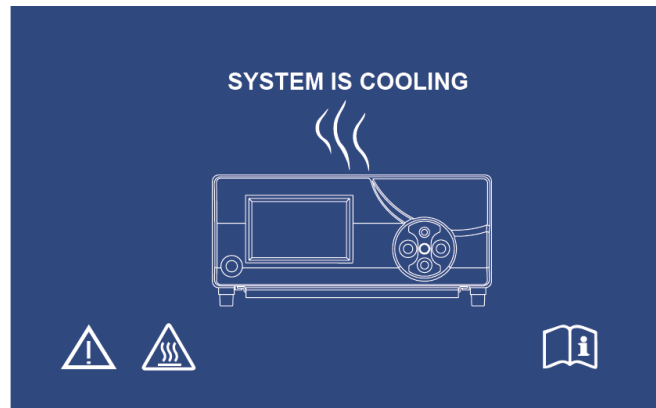
## 6. CLEANING AND DISINFECTION

Both the camera head and camera console may be cleaned, but not sterilized. Remove from power source before any cleaning is performed.

**CAUTION: Never immerse or sterilize the camera console as this will damage the camera console and void the warranty.**

The camera head, camera console and headband can be wiped down with commercially available cleansers commonly used for disinfection of electronic equipment in hospitals such as ethyl or isopropyl alcohols, disinfecting sprays containing quaternary ammonium compounds, or hydrogen peroxide. The lenses should be cleaned only with lens tissue, available in any camera store; follow the directions on the package.

**Do not use strongly caustic or acidic cleansers such as “Clorox” hypochlorite bleach, ammonia, muriatic acid or similar products. Do not use acetone, methyl ethyl ketone or halogenated / chlorinated hydrocarbon solvents or cleansers containing any of these restricted compounds.**

Apply cleaning agents by light spray or dampened towels.

**Do not pour liquids onto the device. Do not allow frank liquids to enter either the camera head or camera console seams or ventilation openings.**

Follow all applicable blood borne pathogen procedures as required by OSHA and/or your hospital, when cleaning and disinfecting the product.

## 7. MAINTENANCE, SERVICING, WARRANTY AND REPAIR

### 7.1. USER SERVICE AND MAINTENANCE

There are no user serviceable parts on the camera head nor the camera console.

Performance of preventative maintenance is not essential. Regular maintenance can, however, contribute to identifying potential problems before they become serious; thus, enhancing the instrument's reliability and extending its useful operating life. Maintenance services can be obtained from your local representative or from the manufacturer.

Defective items or equipment are to be serviced and repaired exclusively by persons authorized by the manufacturer. All repair work shall employ original manufacturer's parts only.

### 7.2. LIMITED WARRANTY

The camera head and camera console carry a 3-year warranty from the date of shipment on workmanship and all defects of material, excluding, fiberoptic cables. Should your product prove to have such defects within three years of shipment, **Sunoptic Technologies**® will repair or replace the product or component part without charge. Should your product(s) need servicing under this warranty, please contact **Sunoptic Technologies**® for return authorization documentation. Please carefully pack the unit in a sturdy carton and ship it to the factory. Include a note describing the defects, your name, telephone number and a return address. Warranty does not cover equipment subject to misuse, accidental damage, normal wear and tear or if transferred to a new owner without authorization from **Sunoptic Technologies**®. This warranty gives you specific legal rights and you may also have other rights that vary from state to state.

### 7.3. REPAIR

You may return your product(s) for repair, shipping prepaid to the factory. Your product will be inspected, and an estimate of repair charges will be submitted to you for approval.

In the US call: (877) 677-2832 (toll free)  
FAX number: (904) 733-4832

## 8. END OF PRODUCT LIFE

This product contains electrical waste or electronic equipment. It must not be disposed of as unsorted municipal waste and must be collected separately in accordance with applicable national or institutional related policies relating to obsolete electronic equipment. The HD camera system must be disposed of according to local laws and hospital practices.

Warning: The camera console contains a lithium coin cell battery. This battery must be disposed of properly according to local regulations.

## 9. SYMBOLOGY

	Manufacturer
	Date of manufacture (YYYY-MM-DD)
	Symbol for "Authorized Representative in the European Community."
	Caution, consult accompanying documents
	Caution: Federal law restricts this device to sale by or on the order of a (licensed healthcare practitioner)
	Consult Instructions for Use
	CE mark
	Not for disposal in general waste
	Caution: Hot Surface
	Caution: Dangerous Voltage
	Product Safety Mark
	Do not use if package is damaged.
	AC Current
	Equipotentiality
	Power on
	Power off
	Type BF
	Protective Earth (Ground)

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